


# Model 8046 Acoustical Enclosure

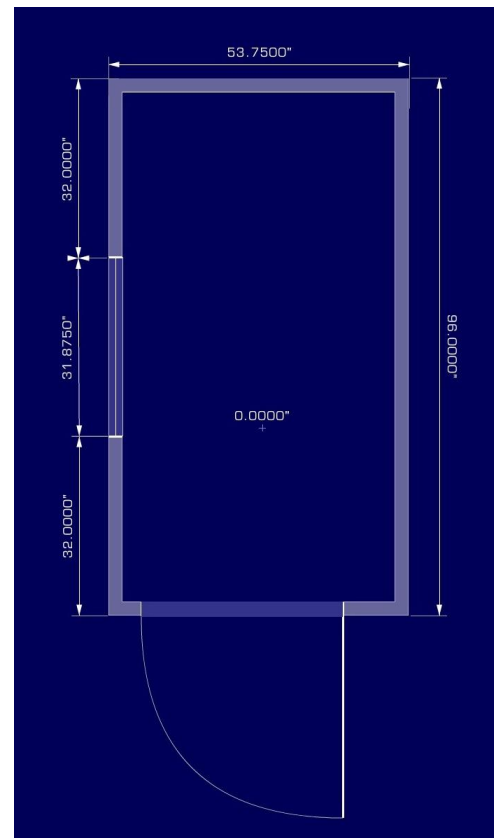

Acoustical Enclosures designed and manufactured in Rock Spring, Georgia.

- Save by buying directly from the manufacturer!
- 2-3/4" thick panels designed to provide the same attenuation as 4" panels allows more interior space than a similar size booth. 18ga solid metal exterior, 18ga perforated metal interior
- Double glazed 1/4" laminated safety glass Window (24"x32") allows live voice testing at a normal voice level
- 36" door opening, Right or Left hinged Door with magnetic seal and compression seal
- 2-1/2" Low Profile Floor with QuietZone Acoustical Floor Mat Base and wall to wall carpet
- 2 LED Light fixtures (low heat, no noise)
- 8 1/4" stereo jacks mounted directly in outer layer of window panel reduces noise transmission through wall and custom placement of cable pass-through, usually at the bottom of the window panel
- Silenced Ventilation with quiet exhaust fan keeps the interior comfortable
- Comparable noise reduction to the major manufacturers and superior to booths constructed of MDF board covered with carpet
- Ready to plug in to a standard outlet after installation

**\$10,900.00** (DELIVERY AND  
INSTALLATION NOT INCLUDED, CALL FOR  
QUOTE)



Telemetrics Corporation, PO Box 129, Rock Spring, Georgia 30739

Telephone: (800)523-8583, Fax: (706)375-7011, e-mail: [bruce@telemetricscorp.com](mailto:bruce@telemetricscorp.com) , [www.telemetricscorp.com](http://www.telemetricscorp.com)